

ARCADIS Avantic

SP

Disposal Instructions

System

Disposal Instructions

© Siemens AG 2005

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Document revision level

The document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

Disclaimer

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

Table of Contents

3

1

Environmental Protection Measures

4

Disposal list

4

2

Changes to Previous Version

7

Disposal list

- The items in the following table are to be disposed of by a suitable disposal company.

Component/assembly	Position		Hazardous material/risk	Quantity
X-ray tube assembly	C-arm, bottom	(1/Fig. 1 / p. 4)	Lead oil implosion of the X-ray tube	1
Tube assembly collimator	C-arm, bottom of the X-ray tube assembly	(2/Fig. 1 / p. 4)	Lead	1
33 cm X-ray image intensifier	I.I. housing, inserted in housing	(3/Fig. 1 / p. 4)	Lead implosion of the image intensifier tube	1
33 cm X-ray image intensifier	C-arm, top lead	(4/Fig. 1 / p. 4)	Lead implosion of the image intensifier tube	1
TFT monitor(s)	Monitor cart, top	(1/Fig. 2 / p. 5)	n.a.	1 / 2
PC	PC in the monitor cart	(1/Fig. 3 / p. 5) (2/Fig. 2 / p. 5)	Lithium	2
Printer	In the monitor cart	(3/Fig. 2 / p. 5)	n.a.	1
Battery	PC	(1/Fig. 3 / p. 5)	Lithium	1
Battery	UPS	(4/Fig. 2 / p. 5)	Lead	1

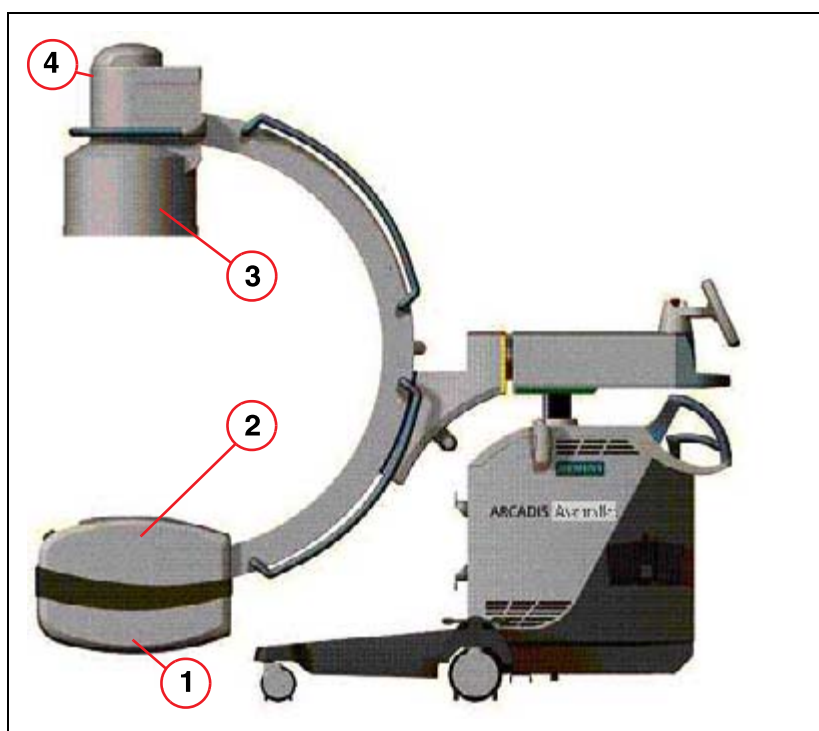


Fig. 1: System overview_ARCADIS Avantic

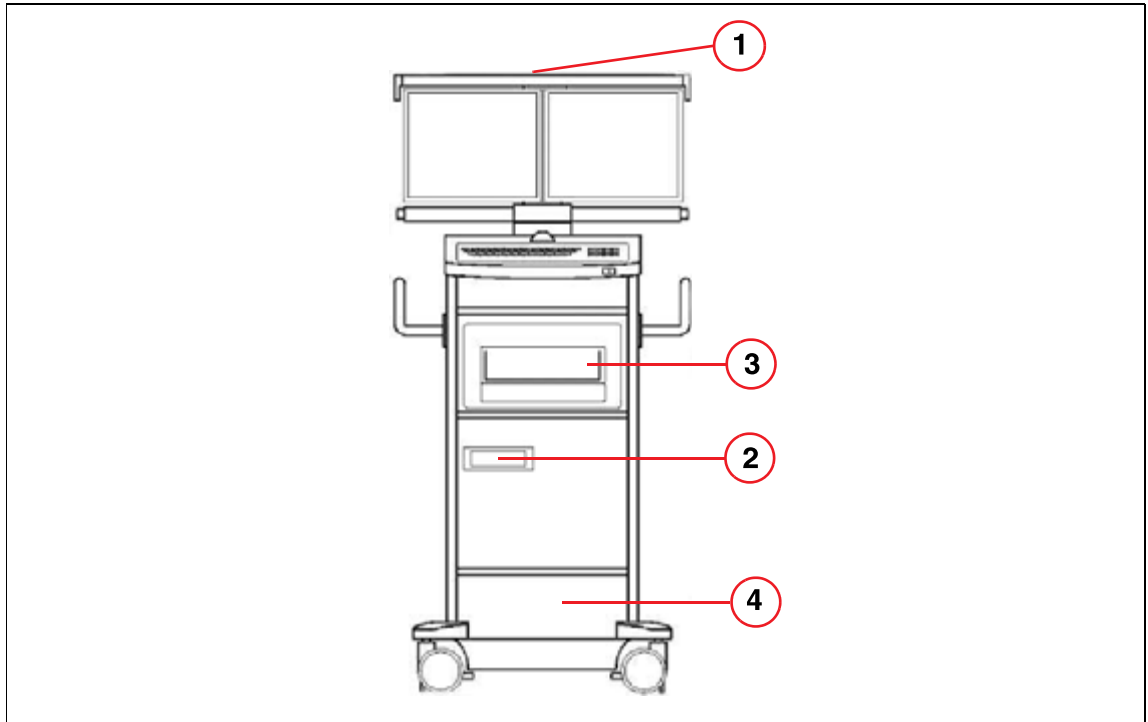


Fig. 2: Monitor trolley

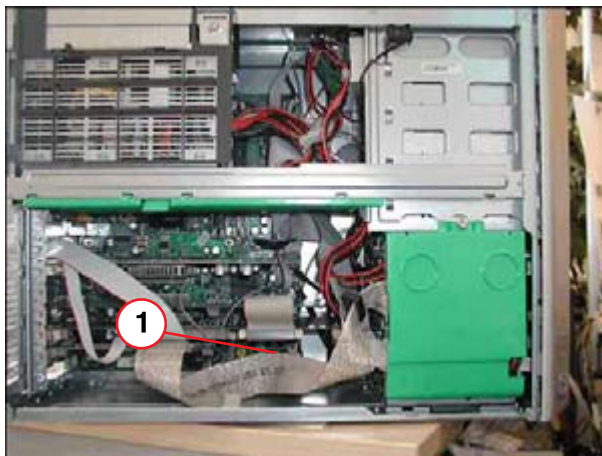


Fig. 3: PC overview
Pos. Lithium battery

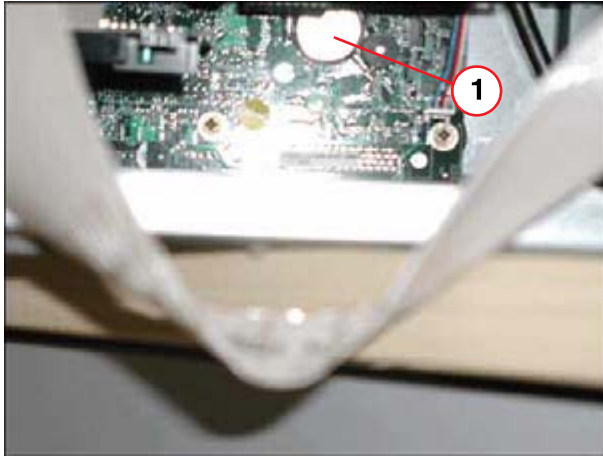


Fig. 4: View of battery

Pos. Lithium battery

n.a., first issue

